

OXYGENIX 2

Wellness Type Hyperbaric Chamber



2.0 / 2.4

ATA working pressure

Sitting Type

Chamber with luxury recline seats

Aluminum

Body, separable modular design

BIBS

100% oxygen breathing system



Spacious interior with luxury leather upholstery, comfortable recline seats, TV and inner control panel.

Spacious & Comfortable Interior


- Wide interior volume
- Comfortable luxury recline seats
- Luxury leather upholstery
- Inner and outer control panels
- Intercom and TV



Specially designed BIBS regulator and masks for comfortable breathing of 100% oxygen.

BIBS Oxygen Breathing System

- Specially designed BIBS regulator
- Comfortable masks for breathing 100% oxygen
- Oxygen and air piping, valves

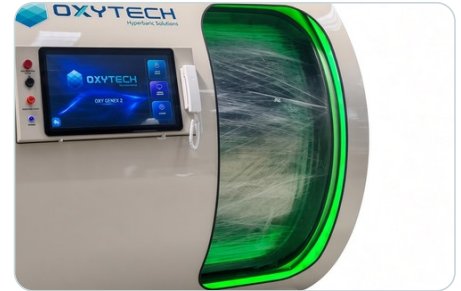

Wide entry door made of acrylic that slides left and right for easy access.

Wide Acrylic Sliding Door

- Acrylic material, wide entry
- Slides smoothly to the left and right
- Acrylic door sliding and locking system

Body & LED Color Options



White body, white LED and entry step



Black body with green LED illumination



Black body, cyan LED, orange seats

TECHNICAL SPECIFICATIONS

Features, Equipment & Dimensions

Chamber Features

- Sitting type chamber
- 2.0 or 2.4 ATA working pressure
- Aluminum body
- Separable modular design
- BIBS oxygen breathing system
- Luxury recline seats
- Intercom and TV
- Oxygen and air piping, valves
- Acrylic door sliding and locking system
- LED illumination
- Interior leather upholstery
- Skirt molding package
- Outer cover
- Step for easy entry into the chamber


Certificates

- Weld certificates
- Non-destructive tests
- Hydrostatic test
- Pressure vessel certificate (2014/68/EU or ASME Section 8, Division 1)

Supplied Equipment

- 2 x Oil-free compressors, 350 LPM
- 2 x Oxygen concentrators, 20 LPM; 93% $\pm$ 3% O₂ purity
- 1 x Cooling system (chillers & radiators)
- Automation software and hardware: special software; outer and inner control panels & screens; pressure, temperature, humidity and O₂ sensors

Dimensions & Weight

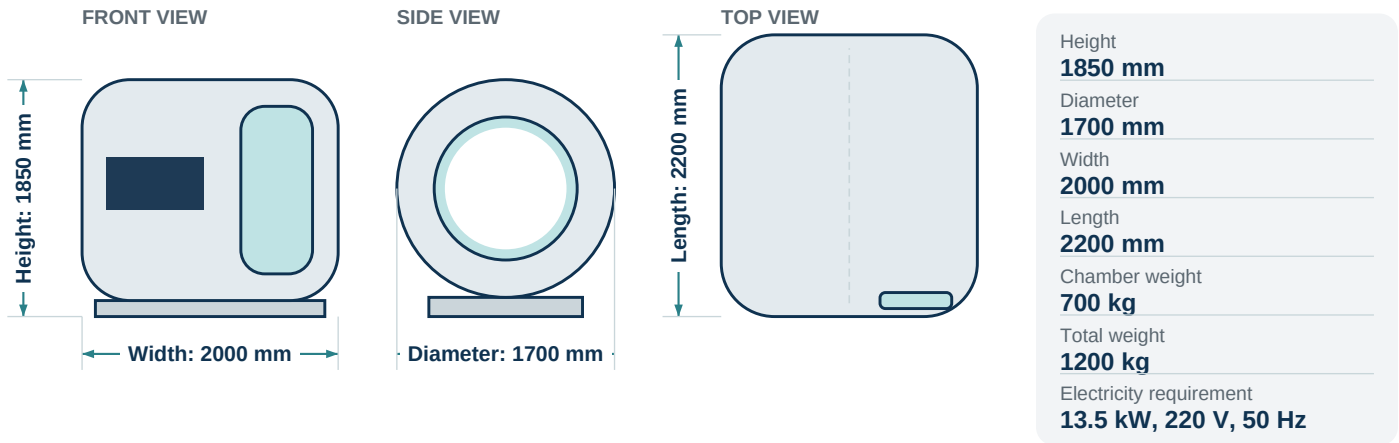

Packaging Details

Wooden Case 1

Dimensions
230 x 235 x 220 cm

Gross weight
1530 kg

Wooden Case 2

Dimensions
135 x 155 x 140 cm

Gross weight
520 kg

CONTACT

Mayer Metal Makina Company

ADDRESS Karadenizliler Mh. Akkol Sk. No: 24/1
Başıskele / Kocaeli / Turkey

EMAIL info@oxytech.com.trWEB www.oxytech.com.tr